**INSTALLATION INSTRUCTIONS** Swiss DolorClast® Cart

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

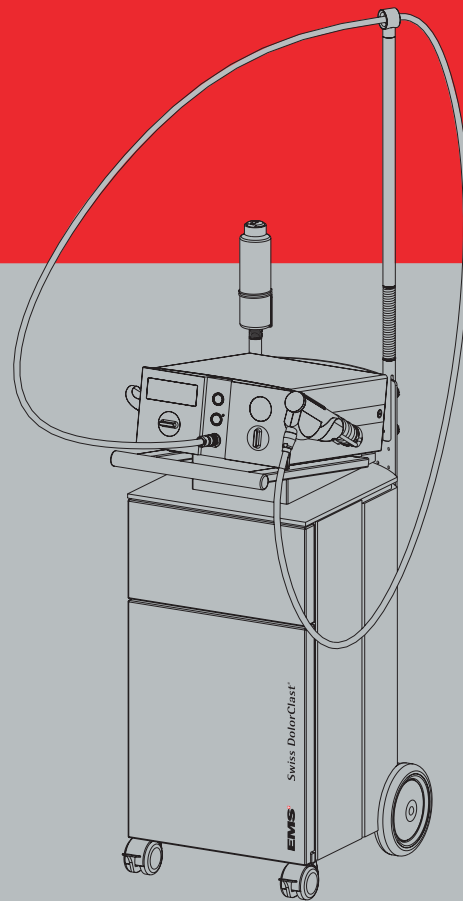
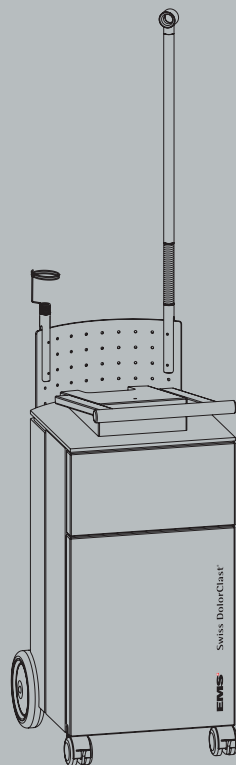
**EMS**   
ELECTRO MEDICAL SYSTEMS
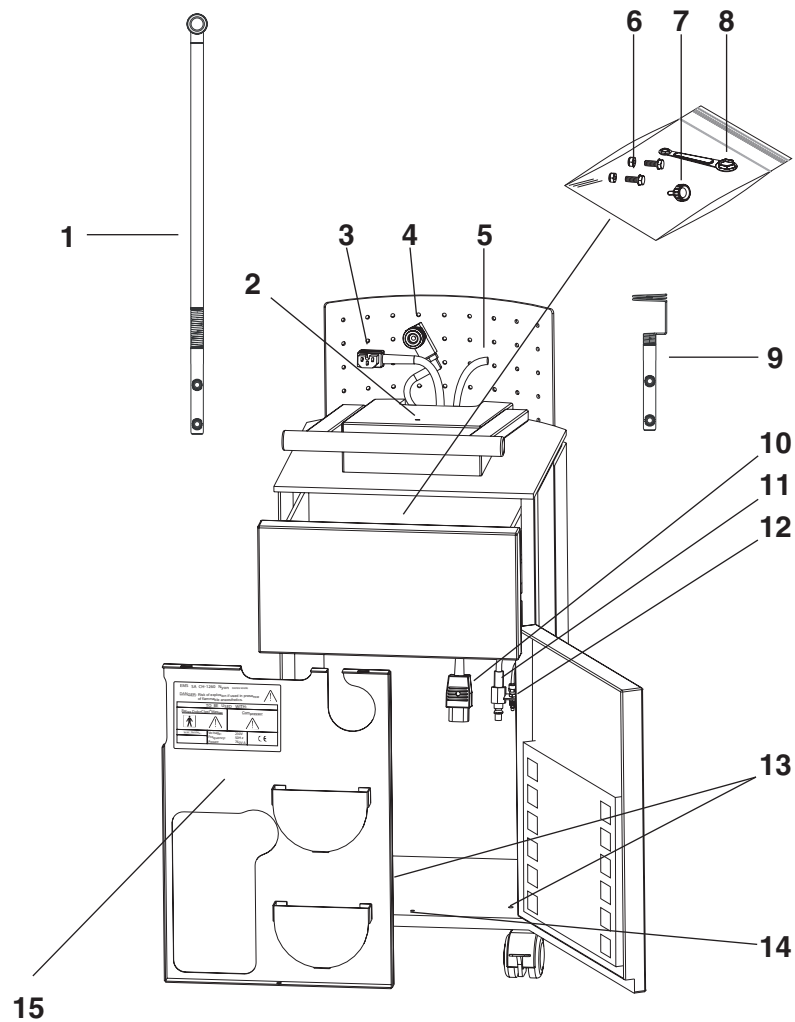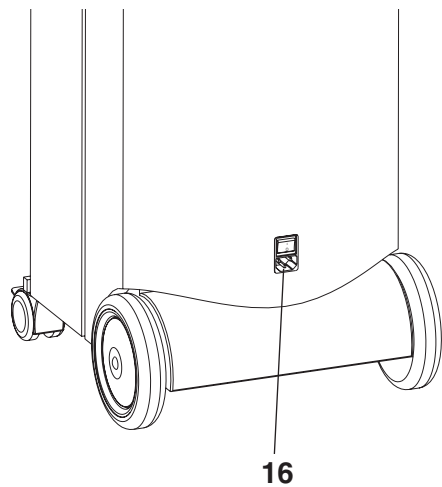
## COMPONENTS

- 1 Carry-cable arm
- 2 Attachment hole - unit
- 3 Power plug - unit
- 4 Air connection - unit
- 5 Connection hose - pneumatic foot pedal
- 6 2 screws + 2 nuts
- 7 Milled screw
- 8 Fitting key
- 9 Bottle holder
- 10 Power plug - compressor
- 11 Connector - air supply - compressor
- 12 Connector - pneumatic foot pedal
- 13 Attachment holes - compressor
- 14 Attachment hole - cover plate
- 15 Cover plate
- 16 Power switch

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## DEAR CUSTOMERS,

Thank you for buying this new EMS product. It meets the highest quality and safety standards.

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### Technical modifications

EMS reserves the right to modify the technique, accessories, operating instructions or contents of the set due to technical or scientific improvements.

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## PLEASE READ THIS BEFORE BEGINNING WORK!

This operating instructions are to ensure the correct installation and use of this product.

Please read these operating instructions carefully as they explain all the most important details and procedures. Please pay special attention to the safety precautions.

Always keep these instructions close at hand.

To prevent injury to people and damage to property, please heed the relevant warnings and remarks. They are marked as follows:

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We would be pleased to answer your questions or listen to your suggestions. We do, of course, provide support in case of technical problems. Please call our EMS service staff or your dealer. The telephone numbers and addresses of the EMS service staff covering your area are listed on the back of these operating instructions.

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EMS

## CONTENTS

Safety precautions.....	2
Installation and setting up.....	6
Preparation of the compressor .....	6
Placing the compressor .....	10
Installation of the unit.....	16
Installation of the bottle holder.....	20
Installation of the carry-cable arm .....	22
Connecting the handpiece cord and placing the gel bottle.....	24
Operating the cart.....	26
Switching On / Off.....	26
General hints .....	28
Cleaning, disinfecting, rinsing, and sterilizing.....	32
Cleaning the cart .....	32
Regular cleaning and maintenance .....	34
Storing the product when out of use for extended periods .....	34
Disposing of the product, accessories and packaging .....	36
Warranty .....	36
Symbols.....	38
Technical data.....	40

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## SAFETY PRECAUTIONS

EMS and the distributor of this product accept no liability for direct or consequential injury or damage resulting from improper use, arising in particular through non-observance of the operating instructions, or improper preparation and maintenance.

 For further details on the preparing, operating and maintenance of this product, please refer to the relevant operating instructions of the unit and handpiece kit.

 Use for the intended purpose only: Before using the product, make sure that you have studied and understood the operating instructions. This also applies to any equipment used with this product. Failure to observe the operating instructions may result in the patient or user suffering serious injury or the product being damaged, possibly beyond repair.

 This Cart is only for use with the Swiss DolorClast® Master (Medical device in compliance with the European Directive 93/42/EEC) and the dedicated compressor. This system is in compliance with the harmonized standard EN 60601-1-1.

Do not connect any multiple plugs or any devices other than the Swiss DolorClast® Master and the compressor.

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Risk of explosion: Do not use this product in the presence of inflammable anaesthetics or gases.

This product must be used only by trained and qualified personnel.

Always examine the product for damage before commencing treatment. Damaged accessories or a damaged unit must not be used and must be replaced. Use original EMS spare parts and accessories only.

Have this product exclusively repaired by EMS or by an authorized EMS partner.

Before each use, always clean, disinfect, rinse and sterilize the different pieces and accessories of this product. Please refer to the information provided in these operating instructions. Non-sterile pieces and accessories may cause bacterial or viral infections.

The cart is unstable, if the compressor is not installed.

This product has been investigated with regard to safety from electrical shock and fire hazard. CSA (Canadian Standard Association) has not investigated other physiological effects.

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
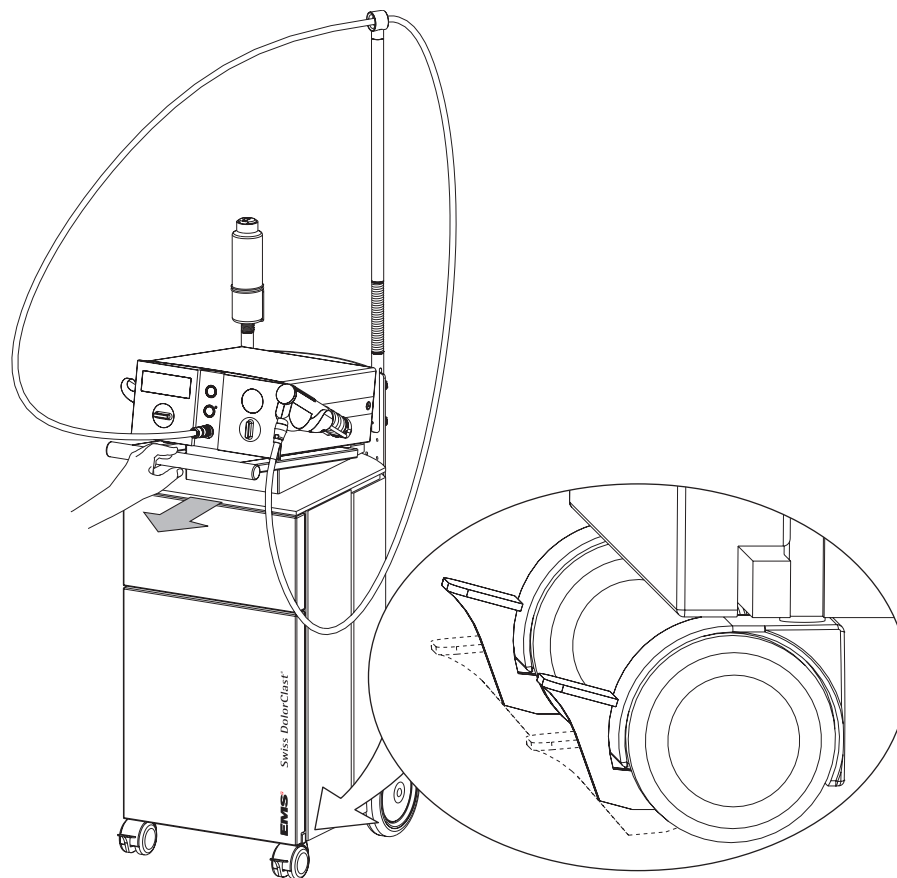
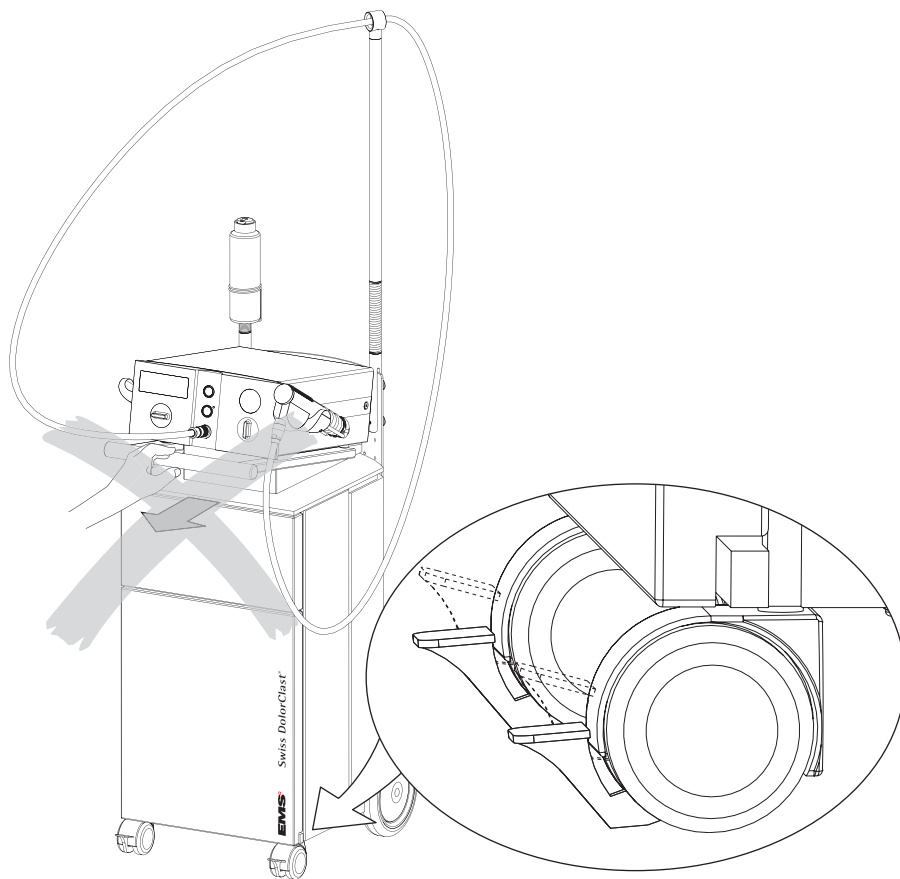
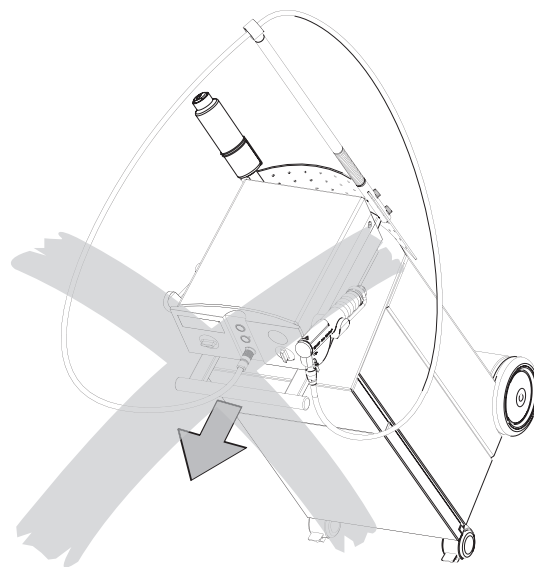
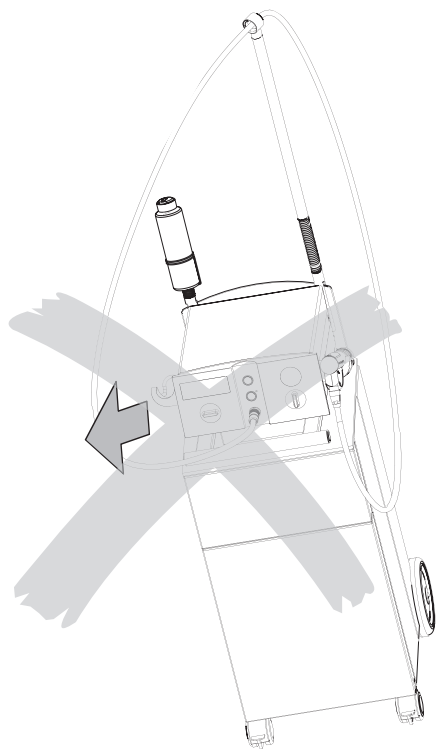
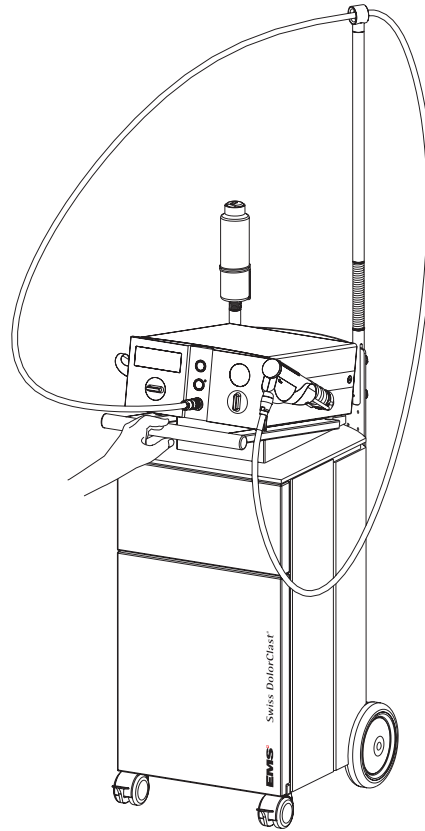
## CLEANING, DISINFECTING, RINSING, AND STERILIZING

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### Cleaning the cart

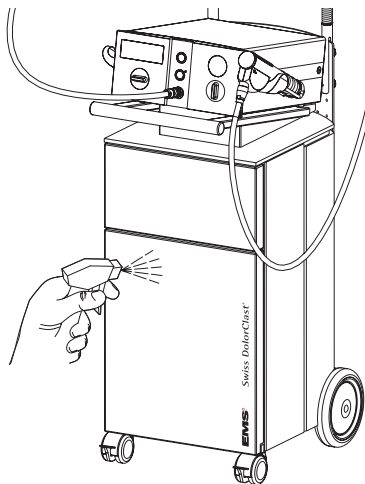
⚠ Clean the cart only with an alcohol-based, commercially available (ethanol, isopropanol), colourless disinfectant. The use of scouring powder or an abrasive sponge will damage its surface.

⚠ The cart is not protected against water spray.

⚠ For detailed information about the cleaning procedure of the unit and the cleaning, disinfecting, rinsing, and sterilizing procedure of the handpiece, please refer to the respective operating instructions.

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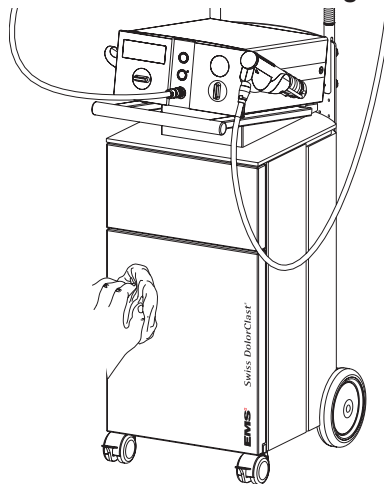## REGULAR CLEANING AND MAINTENANCE

For detailed information concerning regular cleaning and maintenance procedures of the compressor, the unit and the handpiece, please refer to the respective operating instructions.

## STORING THE PRODUCT WHEN OUT OF USE FOR EXTENDED PERIODS

Keep the original packaging until the product is to be disposed of permanently. You can use it for shipping or storing your unit at any time.

Should you wish to put your product out of use for an extended period of time:

- Proceed as described in the chapter "Cleaning, disinfecting, rinsing and sterilizing"
- Pack the product and all accessories in the original packaging

Please refer to the storage and transport conditions in the "Technical data".

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## DISPOSING OF THE PRODUCT, ACCESSORIES AND PACKAGING

The product, its accessories and its packaging do not contain any substances that are dangerous for the environment.

Should you wish to definitely dispose the product, please comply with the regulations which apply in your country.

 For further details concerning the disposing of the compressor, please refer to the operating instruction.

## WARRANTY

The warranty will be valid for one year from the date of purchase of your product and accessories.

Damages due to non-adherence to the operating instructions or wearout of parts are excluded from warranty.

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EMS MEDICAL GmbH, Obere Laube 44, D - 78462 Konstanz
EMS Corp. USA